

Tasso +
Instructions For Use
AGE 18 AND OLDER

tasso+



Device Description

Tasso+ is a sterile, single-use blood collection device. It adheres to the skin with a mild adhesive and, when activated, a lancet creates a small puncture while a gentle vacuum draws capillary blood into a collection tube. Once the tube is full or the collection time ends, the device is removed, and the tube is capped for transport. The lancet retracts automatically for safe disposal.

Precautions

- For use only on a single patient.
- Discard the entire device after use
- Do not use after expiration date
- Single use only.
- Keep out of reach of children.
- Prescription use only.
- Do not resterilize.
- For external use only.
- Do not use if device packaging has been opened or damaged.
- Use while seated as fainting may occur with any blood sampling procedure.
- Minor bruising, residual marks, or scarring may occur at the sample collection site.
- Patients taking blood thinners may experience prolonged bleeding.

Intended Purpose

The Tasso+ is a single-use blood lancing device intended for obtaining microliter capillary whole blood samples.

First time using Tasso?



Watch the instructional video:
www.tassoinc.com/tassoplus-video

WARNINGS: Not intended for more than one use. Do not use on more than one patient. Improper use of blood lancets can increase the risk of inadvertent transmission of bloodborne pathogens.

1 Prepare

Step 1

Wash hands before use. Work on a clean, flat surface with space to lay out the kit.



Step 2

Device will be placed on upper arm - for best results remove any excess hair.



Step 3

Take out all items from kit. Remove device from tray. Unwrap bandage and gauze and set aside.




Step 4

Twist off the cap from the tube. Set aside DOB sticker label for later use.



Step 5

Press tube into device until snug. Be mindful of the sticker and position the tube so you can watch it fill.

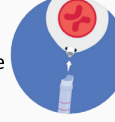

Step 6

Clean area with alcohol pad and allow to dry. Use the alcohol pad to clean the heat pack.



Step 7

Follow heat pack instructions to activate. If pack isn't available, use other warming method such as hand.



Step 8

Rub heat pack or hand up and down shoulder for 4 minutes to warm and increase blood flow



Step 9

Remove clear cover from red button.



2 Collect

Step 10

Peel tab away from device and stick device to upper arm.



Hang arm straight down at side.

Remain seated during collection to reduce the risk of dizziness or fainting.

Step 11

Press button all the way down ONCE and release. Set a timer for 5 minutes. Watch tube fill using mirror as needed.



Push firmly - some resistance is normal and you may hear a faint click.

Step 12

Blood may not appear for first 1-2 minutes.



Remove device at 5 min total or as soon as the blood reaches top fill line.

3 Finishing Collection

Step 13

Have gauze and bandage ready. Peel the device off from one side.

Step 14

Press gauze on site and apply bandage vertically with light pressure.



Step 15

Remove tube from device with a slight twist and pull down.



Step 16

Snap cap fully onto tube until it clicks. Ensure cap is securely on tube.



Gently Invert tube 8-10 times. **Do NOT** shake tube.



4 Labeling & Transport

Step 20

Write your Date of Birth on the sample tube label.

DATE OF BIRTH:
MM / DD / YY.

Date of Birth
MM / DD / YYYY



Please note that failure to do this will result in your sample being rejected at the lab.

Step 21

Complete the included form with your **Full Name, Date of Collection, DOB, and signature.** Then follow the provided shipping instructions to mail your sample.



Failure to complete the form will result in sample rejection.